

BAX[®] System X5

PCR Assay for *E. coli* O157:H7

The BAX[®] System X5 offers the same fast, accurate, easy-to-use pathogen detection solution that customers have come to expect from the BAX[®] System method, but in a smaller, lightweight construction. This PCR assay detects *E. coli* O157:H7 from a wide variety of enriched samples (meats, vegetables, fruits and raw milk), using automated sample preparation in the Hygiena[™] Thermal Block and automated amplification and detection in the BAX[®] System X5 instrument.



QUA 18/04 - 03/08
Alternative Analytical
Methods for Agribusiness
<http://nf-validation.afnor.org/en>

Features & Benefits:

- Clear yes-or-no results in as little as 12 hours for select food samples
- Compatible with other BAX[®] System X5 assays for efficient processing
- Carefully designed primers target specific genetic sequences possessed only by the target organisms
- Minimal components and simplified workflows to maximize efficiency and ease-of-use
- Validated to perform as well or better than standard reference methods
- Internal controls included to validate results even in absence of target
- Flexible protocols available to meet your unique workflows

Validations, Certifications and Adoptions:

- **AOAC Research Institute**
Performance Tested MethodSM #050501
Validated on raw ground beef, beef trim, spinach and lettuce
- **NF VALIDATION certificate granted by AFNOR Certification**
QUA-18/04 - 03/08 - Certified according to AFNOR validation rules for raw beef meats, raw pork, ovine and chicken meats, ready to eat and ready to reheat meals, raw milk and fruits and vegetables

keydiagnostics


T: 02 8212 4074 F: 02 9423 6992
info@keydiagnostics.com.au
www.keydiagnostics.com.au
PO Box 1038, Gymea, NSW, 2227

Hygiena Product Code	Legacy Order Code	Description	Quantity
KIT2022	D15407214	BAX [®] System X5 PCR Assay for <i>E. coli</i> O157:H7	64 per kit




Find support documents, instructional videos, and more at www.hygiena.com


BAX® System Protocol

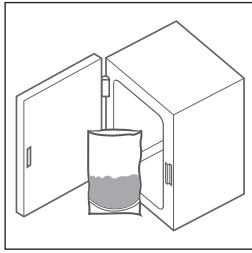


Enrich Samples.

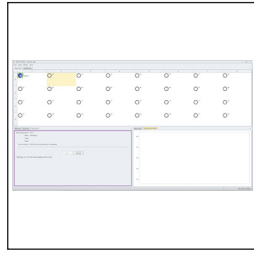


Create rack file and warm up cyclers.

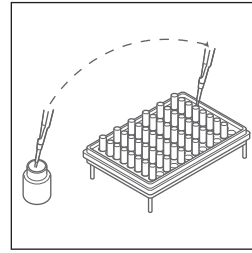


Mix protease with lysis buffer and transfer 200 µL of lysis reagent to cluster tubes.

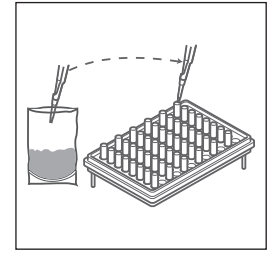


Transfer 20 µL sample enrichment to cluster tubes.

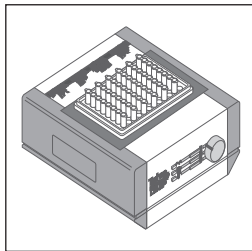

Place samples on automated thermal block for lysis and cooling.

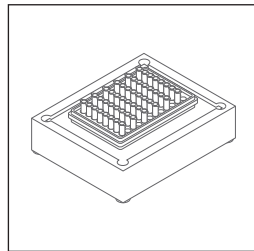

Transfer 50 µL of lysed sample to PCR tubes in cooling block.

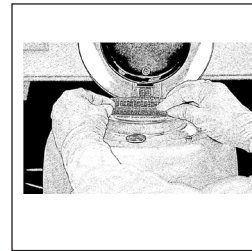

Place sealed PCR tubes in cyclers and immediately click "NEXT" to run program.



Review MP Express results. (Add 60 minutes for standard MP results.)

Related Products

BAX® System MP Media

Available enrichment media for customers looking to take full advantage of the rapid time-to-result and ease-of-use offered by select BAX® System *E. coli* and *Salmonella* assays.

StatMedia™ Soluble Packets

Gamma-irradiated BAX® System MP Media in convenient, water-soluble packets for reduced mess and preparation. Simply drop in pre-warmed sterile water and mix with sample.

BAX® System X5 PCR Assay for *Salmonella*

Detects *Salmonella* species from standard enrichments of a variety of food and environmental samples with the same accurate and reliable DNA-based pathogen detection as the BAX® System Q7 with a smaller footprint.

BAX® System X5 PCR Assay for *L. monocytogenes*

Detects *Listeria monocytogenes* from a wide variety of enriched samples (meats, seafood products, dairy products, fruits and vegetables) with the same accurate and reliable DNA-based pathogen detection as the BAX® System Q7 with a smaller footprint.

Hygiena Product Code	Legacy Order Code	Description	Quantity
MED2003	D12404925	BAX® System MP Media	2.5 kg tub
MED2016	D12745725	StatMedia™ Soluble Packets	20x5x33.75g
KIT2025	D15407187	BAX® System X5 PCR Assay for <i>Salmonella</i>	64 per kit
KIT2023	D15407201	BAX® System X5 PCR Assay for <i>L. monocytogenes</i>	64 per kit



Find support documents, instructional videos, and more at www.hygiena.com/bax