

# GLUTENTOX STICKS

## GlutenTox Stick Kit Difference



## Rapid Gluten Detection Kit

### The GlutenTox Difference

GlutenTox Sticks is a rapid and sensitive lateral flow test for the detection of immunotoxic gluten in raw ingredients, finished products and working surfaces.

At a detection 3ppm, it is one of the most sensitive kits on the market, so you can be well assured of your gluten free claims.

### GlutenTox Advantages

- Semi Quantitative results
- Specific to the toxic fraction of gluten
- Detection limit 3ppm
- Surface results in 15 minutes
- Product results in 1 hour
- Suitable for hydrolysed foods
- No false positives or negatives

The GlutenTox kit utilises a new anti-Gliadin antibody (G12) for the detection of the 33 mer peptide toxic fragment of gluten, the primary initiator of the inflammatory response to gluten in Celiac Disease.

#### SENSITIVITY

Detection limit 3ppm (15ng/mL) after 10 mins


#### SPECIFICITY

Specific to prolamins of,

- Wheat (Gliadin)
- Barley (hordein)
- Rye (secalin)
- Oat (avenin)

Doesn't cross react with,

- Soy
- Rice
- Corn

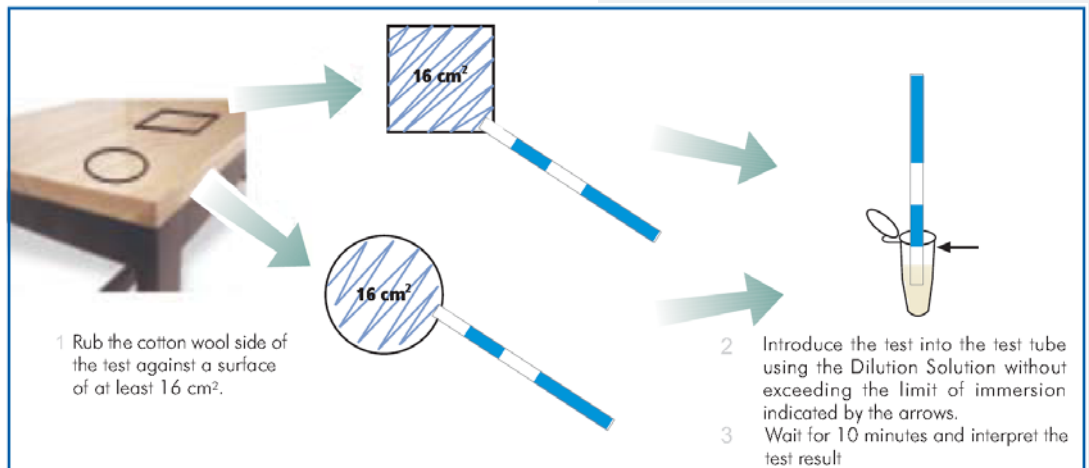

Tel : 02 8212 4074 Fax : 02 9423 6992  
PO Box 1038 Gymea NSW 2227  
Email: info@keydiagnostics.com.au

[www.keydiagnostics.com.au](http://www.keydiagnostics.com.au)

