

Introducing Puritan's Environmental Sampling line

The result of nearly a century of experience –
a connection to your markets – an innovative
R&D team – and super-stringent QC.



EnviroMax[®]
(dry)

EnviroMax Plus[®]
(pre-moistened)

Puritan[®]

Quality Medical Products Since 1919

KEY
diagnostics



Tel : 02 9427 7686

Fax : 02 9423 6992

www.keydiagnostics.com.au

PO Box 1038, Gymea NSW 2227

100% US-MANUFACTURED

EnviroMax® and EnviroMax Plus® – the ultimate in environmental sampling from the swab experts at Puritan.

We spent time **in the marketplace** with different end users, from **food production** to **cosmetics** to **pharmaceutical firms**. We saw how you work – listened to what you need – then went to our R&D team to design best-in-class environmental swabs.



The EnviroMax® and EnviroMax Plus® are oversized, perfect for sampling surfaces or large equipment parts. Count on Puritan for consistent quality, time after time.

1

Quick-turn-open, leak-resistant cap. No worry about media solution leaks or loss of post-collection contents.

2

Swab securely attached to cap. Gives you an excellent grasp for maneuvering the swab while testing.

3

Rigid paddle inside foam. Ideal for reaching into crevices with precision and control.

4

Available pre-moistened or dry. Your choice, depending upon your application.

5

Sterile wrapped. Features easy edge-peel wrapper.

6

White handle and swab tip. Makes it easy to view any post-collection colorimetric changes.


7

Maximized collection and elution. This was the goal of our tip materials and shapes – to give you the best environmental swabs.

8

Swab tip shape options. We know that testing surfaces vary – so we offer tip options and do custom orders!


Products on this page shown at 70% of actual size.

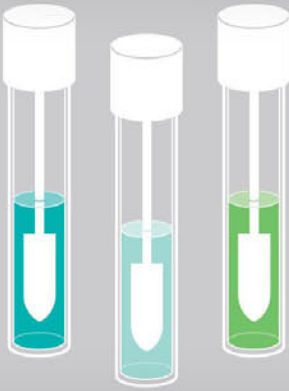

WHITE HANDLE & TIP, CLEAR VIAL

There are often times when you do a test to get an immediate result – or at least an indicator that a particular sample should go to the lab for more extensive testing.


That's why all the swabs in our environmental line have white handles and tips and are packaged in a clear plastic vial. Reactions – even small, early reactions – are visible, saving you time and effort. Count on us for reliable, repeatable results, clearly seen.

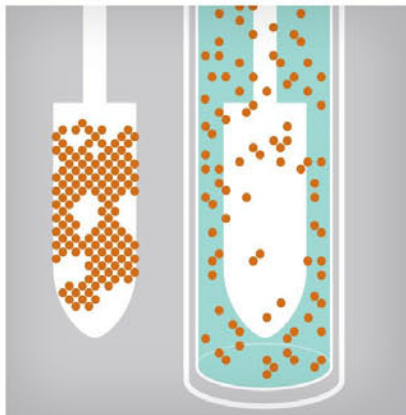

CATCH & RELEASE

It's just common sense. The more sample you can pick up the more accurate the results. That's why Puritan worked so hard to maximize the collection surface for our environmental swabs. The tip fibers act like a "fishing net" to catch and hold more than enough material for accurate testing.

At the same time, these fibers readily release sample material once the swab is inserted into media. Catch and release – what Puritan swabs are known for.

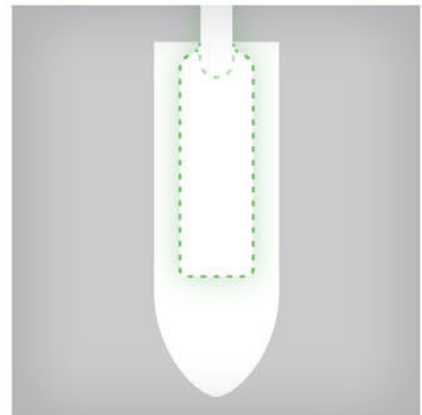

PADDLE STRUCTURE

The EnviroMax® and EnviroMax Plus® have oversized tips with a rigid paddle-like structure within the foam.

This gives the tip stability and makes it easy to reach into crevices where just foam may bunch up, crush, or fold over.

These were inspired by a customer in the food production industry where large pieces of equipment with varying surface terrain need sampling.

Why Puritan environmental swabs?

They are **easy to use** – ideal in almost any environment. You can choose options that give you **instant results**, perfect for spot checks of virtually any surface.

As with all Puritan products, you can count on the **highest quality**, with **rigorous QC** and oversight, to give you **peace of mind** knowing that your test results will be **accurate** and **repeatable**.



While Puritan is relatively new to the environmental sampling

field, we bring almost a century of **experience** to the single-use products we develop for you.

While many things have changed – new technologies, materials, and processes – some things are the same as they were in 1919. We design, manufacture, package, and ship products from our Guilford, Maine, **USA facility**. Our products are readily **available** with no waiting. And our **customer service** is second to none.

ENVIROMAX® (DRY) AND ENVIROMAX PLUS® (PRE-MOISTENED)

This part of our environmental line includes larger swabs – you can specify a dry swab kit or one with a pre-moistened swab. These are an excellent choice for sample collection on larger surfaces. Products below shown at 45% of actual size.

	Item No.	Item Description	Packaging pack, weight	Tube length	Swab length	Tip length x diameter
	non-sterile sterile					
	88010-PF	Pointed, dry foam tip, heat staked to polypropylene handle.	25/1, 2 lbs	152.91 mm 6.02"	147.32 mm 5.8"	56.388 x 15.545 mm 2.22 x 0.612"
	25-88010 PF					
	88010-PFB		10/25/1, 20 lbs			
25-88010 PFB						
	25-88050 PF	Pointed, pre-moistened foam tip, heat staked to polypropylene handle. Individually wrapped.	25/1, 2 lbs	152.91 mm 6.02"	147.32 mm 5.8"	56.388 x 15.545 mm 2.22 x 0.612"
	25-88050 PFB					



Headquarters:

1430 West McCoy Lane
 Santa Maria, CA 93455
 800.266.2222 : phone
 805.346.2760 : fax
 Sales@HardyDiagnostics.com
 www.HardyDiagnostics.com

Distribution Centers:

Santa Maria, California
 Centralia, Washington
 Salt Lake City, Utah
 Phoenix, Arizona
 Dallas, Texas
 Springboro, Ohio

Lake City, Florida
 Albany, New York