

3M Infection Prevention Solutions
3M™ Comply™ Lead Free Indicator Tape for Steam Sterilisation

Simple
Reliable
Lead-Free



The industry standard for exposure monitoring. Now lead-free!

- Shows at a glance whether packs have been exposed to a steam sterilisation process
- Accepts hand-writing and labeling
- Does not require hazardous waste disposal



Seals securely



3M Infection Prevention Solutions

Innovation
On A Mission

New lead-free construction – same trusted adhesive performance

The simple, convenient, sustainable way to check packs for sterilisation exposure

3M™ Comply™ Steam Indicator Tapes are external chemical indicators used to securely seal the wrapping on instrument packs. They provide the steriliser operator with visual assurance that each package has been exposed to the sterilisation process, without having to open the pack.

Available in two adhesive types – one for all wraps, another specifically for disposable wraps – Comply Steam Indicator Tapes feature a stretchable backing that minimises tape “pop-off” during sterilisation. And they can be written on or labeled with pre-printed labels, such as 3M™ Comply™ Sterilisation Load Labels.

Healthcare facilities around the world have relied on these easy-to-use, economical indicator tapes for more than 60 years. Now, thanks to our new lead-free construction, you can enjoy that same level of trusted performance – without the hassles and expense of managing and disposing potentially hazardous waste.




3M™ Comply™ 1322 Steam Indicator Tape is designed for use on all wraps. Provides a secure seal and is easy to remove from reusable wraps.

3M™ Comply™ 1355 Steam Indicator Tape is specifically designed for use with disposable wraps. It features an aggressive adhesive to secure packs during and after processing.



The chemical indicator ink, printed in a diagonal striped pattern, changes color or turns darker when exposed to steam.

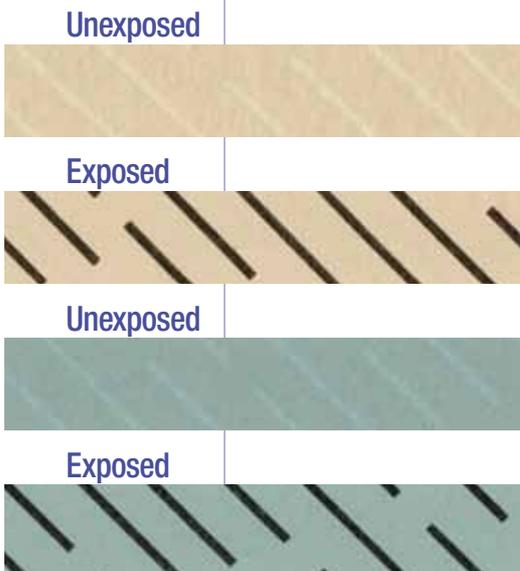



...taking you a step beyond Best Practice

Lead is considered a “hazardous waste” under the Federal Resource, Conservation & Recovery Act (RCRA), which is the basis of U.S. EPA regulations for hazardous waste.

Congress has authorised the states to carry out many of the functions of RCRA through their own hazardous waste programs. For example, nineteen U.S. states now prohibit the sale or distribution of packaging that intentionally contains lead or other toxic metals. The purpose of these laws is to prevent the use of toxic heavy metals that enter landfills, waste incinerators, recycling streams, and ultimately, the environment.

The problem is, hospitals are full of products that could be considered hazardous waste once they are thrown away – and often you might never even know that a product contains lead. For example, lead has been used for decades as an ingredient in the specialised inks used in the manufacture of hospital steriliser indicator/autoclave tapes and cards. When discarded, these can be considered toxic hazardous wastes, depending on the amount and concentration of lead that is disposed of at any given time.

Although lead occurs naturally in the environment, human activities have increased levels of lead in the biosphere more than a thousand-fold in the past 300 years. Much of this increase occurred in the second half of the twentieth century because of the use of lead compounds in petrol, a practice which has ceased in Australia (except under special licence). Lead also enters the environment as a consequence of the mining and refining of lead, industrial and manufacturing activities associated with products which contain lead, and the inappropriate disposal of waste materials containing lead.

As concern about the long-term health effects of lead and other toxic metals continues to grow, it is almost certain that regulations on the use and disposal of lead-containing products will become increasingly restrictive, complex and costly for healthcare facilities to follow. That’s why standardising on lead-free products, such as 3M™ Comply™ Steam Indicator Tapes, make such good sense – by helping to simplify operations, avoid the risk of penalties for non-compliance, and enhancing your reputation for good environmental stewardship.

Reference: Australian Government. National Health and Medical Research Council (NHMRC) Information Paper, August 2009.

“An Information Paper for Practitioners and Policy Makers – Blood lead levels for Australians”

Why it is so important to “get the lead out.”

Lead is a soft, dense, grayish metal that has been used for thousands of years in the manufacture of products ranging from plumbing pipes to cosmetics.

The problem with lead is that, in sufficient quantities, it is poisonous, and can cause a variety of serious health effects, including neurological damage, blindness, kidney failure – and can even lead to death. It is especially dangerous to young children.

Increasingly-stringent regulatory actions, such as the elimination of lead in petrol and household paints, have reduced the amount of lead that reaches the environment; however, lead can still be found in a wide range of products, from fishing tackle and automobile wheel weights to car batteries. If improperly disposed of, lead from these products can leach into the air, water and soil – creating conditions that are potentially lethal to people and aquatic life.

Best Practices

- Use an external chemical indicator on the outside of each package
- If a package allows for visual inspection of an internal indicator (such as those with paper/plastic packaging), an external indicator is not required

Standardise to the Core

3M™ Comply™ Steam Indicator Tapes are one of the “Core Four” 3M sterilisation products. By standardising on the Core Four, and following AAMI and AORN recommended standards and practices, you can be confident that you are doing your best for your facility and its patients.



Ordering Information

Class 1: Conforms to ANSI/AAMII/ISO 11140-1:2005

Catalog Number	Product Name	Size	Rolls/Case
1322-12MM	3M™ Comply™ Lead Free Steam Indicator Tape (designed for use on all wraps)	12mm x 55m	42
1322-18MM		18mm x 55m	28
1322-24MM		24mm x 55m)	20
1355-18MM	3M™ Comply™ Lead Free Indicator Tape (designed for use on disposable wraps)	18mm x 55m	28
1355-24MM		24mm x 55m	20

Note: Contains dry natural rubber.

Indicator Tape Dispenser

Catalog Number	Product Name	Description	Items/Box	Boxes/Case
C22	Heavy Duty Tape Dispenser	Holds 2 rolls	1	2
M52	Heavy Duty Tape Dispenser	With Tabber	1	1



3M Health Care

3M Australia Pty Limited
 ABN 90 000 100 096
 Bldg A
 1 Rivett Road
 North Ryde NSW 2133
 1300 363 878
www.3M.com.au/healthcare

3M New Zealand Limited
 94 Apollo Drive
 Rosedale
 Auckland 0632
 0800 80 81 82
www.3M.com./healthcare

3M is a trademark of 3M Company
 © 2011. All rights reserved
 PB5185/0611