



MiniMix[®] 100 CC[®]

Lab blender

Simple, Accurate, Space-saving

The most compact blender in the world, BagMixer® 100 MiniMix® is now available in Click & Clean® version!

MiniMix® 100 automatically blends any solid sample in a sterile bag and allows safe and easy extraction of bacteria. You can pipet and analyse the filtrate without any risk of alteration or cross contamination. **This is the safest and most economical sample preparation method.**

The new MiniMix® 100 CC® incorporates the Click & Clean® and side-by-side paddle stopping systems, already highly acclaimed on BagMixer® 400 CC® and the adjustable blending strength system .




1 Multifunction screen

For easy control of the blending parameters

Blending time


Blending speed

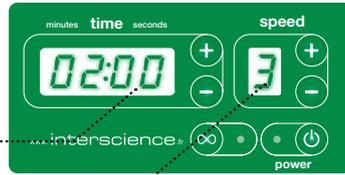

2 Q-Tight® system


A special closing system that avoids every risk of leaking.

Leak-tight

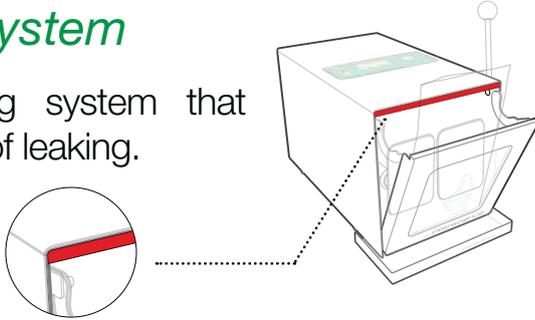

3 Adjustable blending strength **PATENTED**

You can adapt the blending chamber to the size of the sample to get OPTIMAL EXTRACTION!

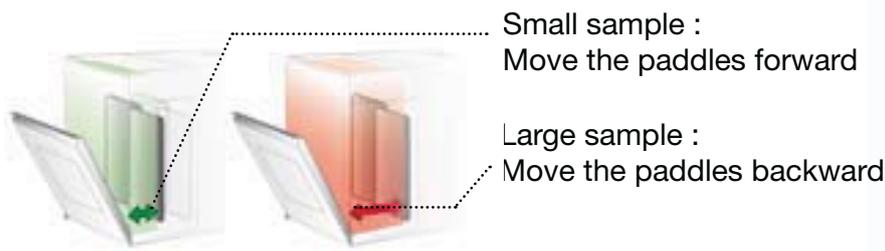


Small sample :
Move the paddles forward

Large sample :
Move the paddles backward


4 Click & Clean® system **PATENTED**

UNIQUE ! TOTAL CLEANING of the blending chamber is now possible !

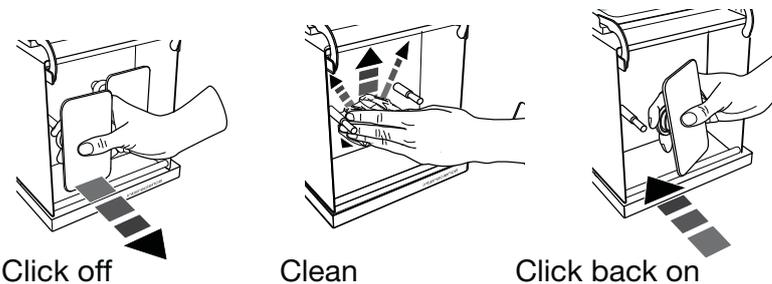

5 Side-by-side paddle stop **PATENTED**

For easy and fast insertion of your sterile bags in the blender.

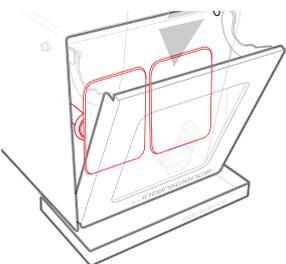

Advantages

- Compact, silent
- User friendly
- Optimal extraction
- Lifetime warranty on window door and shock absorbers
- Stainless steel
- Q-Tight® : Leak-tight
- Easy and total cleaning

Applications

- Food microbiology
- Pharmaceutical
- Veterinary, clinical (biopsies, organ sampling...)
- Molecular extraction (proteins, lipids, DNA...)
- Toxin research (HPLC, GCMS...)

Specifications

- Useful capacity : 5 - 100 ml
- For samples : 0.5 g - 10 g
- Adjustable time : 1 sec. à 59 min. countdown & open running (∞)
- Speed calibrated by microprocessor (3 - 12 strokes/sec.)
- Side by side paddle stop
- Q-Tight® seal : thickness of 12mm

- Stainless steel
- Drip tray
- Adjustable blending strength
- Click & Clean® system : removable & autoclavable paddles
- Multifunction screen
- Lifetime warranty on window door & shock absorbers

- Pwr. : 220 / 100 - 110 V – 50 / 60 Hz
- Dim. (Lxlxh) : 37 x 20 x 20 cm
- Net weight : 13 Kg

Ref : 013 230



WEEE Compliant
Directive 2002/96/EC

RoHS Compliant
Directive 2002/95/EC



Please return this form to
interscience or to your local
distributor and indicate your interests.



Also available

MiniMix® 100 ml blender

- Quotation
 Demonstration

I would like to receive more info on :

- BagMixer® 400 ml blender
 JumboMix® 3500 ml blender
 BagFilter® filter bags
 BagTools® accessories
 Gravimat® gravimetric dilutor
 Spiral® DS+ plater
 Scan® 100, 500 & 1200 colony counters

Name.....

Company.....

Address.....

Zip/City..... Country.....

Phone..... Fax.....

email.....

web.....



Tel : 02 8212 4074
Fax:02 9423 6992
PO Box 1038
GyMEA NSW 2227
info@keydiagnostics.com.au
www.keydiagnostics.com.au

interscience

Worldwide: 30. Ch Bois Arpents, 78860 St Nom FRANCE
Tel. +33(0)1 34 62 62 61 / Fax. +33(0)1 34 62 43 03
info@interscience.fr / www.interscience.fr

USA/Canada: 301 Winter Street, Hanover, MA 02339, USA
Tel. +1 781 829 9026 / Fax. +1 781 829 9027
info@intersciencelab.com / www.intersciencelab.com



MiniMix® 100 VP

- Adjustable speed (6-9 strokes/sec)
- Adjustable time (10 - 360 sec. or ∞)
- Available in CoolMix® version

Get the best results with...



Gravimat® Gravimetric dilutor

Automatically dilutes with the appropriate mass of solvent.



BagPage® 80 & 100 Full page sterile filter bags

Robust & food compatible, multilayer.



BagLight® 100 Plain sterile bags

Made of PolySilk®, extremely robust. For samples that do not require any filtration.



RollBag® 450 Sampling bag

For collecting, sampling & preservation of any sample. Stands up alone.



BagOpen® 100 Bag opener

Easy opening, extremely stable. Resistant to cleaning agents.



BagRack® 100 Storage rack

Stainless steel, autoclavable.



BagClip® 100 Closing clips

100% watertight & waterproof, closes any kind of bag.